

# EmStat4s™

POTENTIOSTAT | GALVANOSTAT | IMPEDANCE ANALYZER (optional)



## Contents

Desktop performance in the palm of your hand .....	3
Supported Techniques .....	4
Measurement Specifications .....	5
System Specifications .....	6
EmStat4S LR EIS Accuracy Contour Plot .....	8
EmStat4S HR EIS Accuracy Contour Plot .....	8
Standard EmStat4S Kit .....	9
PSTrace: Software for Windows .....	10
EmStat4S works with MethodSCRIPT™ .....	12
Software Development Kits for .NET .....	13
PStouch: App for Android .....	14

➤ See for more information:  
[www.palmsens.com/emstat4s](http://www.palmsens.com/emstat4s)



## Desktop performance in the palm of your hand

The EmStat4S is a portable, USB-powered Potentiostat, Galvanostat, and optional Frequency Response Analyser (FRA) for Electrochemical Impedance Spectroscopy (EIS). The EmStat4S Low Range version is great for applications that require measuring low currents down to picoamps, whereas the High Range version is very suitable for applications that need a maximum current of up to 200 mA.



Main differences between the Low Range (LR) and High Range (HR) versions:



### EmStat4S LR

### EmStat4S HR

▪ potential range	±3 V	±6 V
▪ max. compliance	±5 V	±8 V
▪ current ranges	1 nA to 10 mA (8 ranges)	100 nA to 100 mA (7 ranges)
▪ max. current	±30 mA	±200 mA
▪ FRA/EIS (optional)	10 $\mu$ Hz to 200 kHz	
▪ electrodes	WE, RE, CE, and ground, 2 mm banana plugs	WE, RE, CE, S, and ground, 2 mm banana plugs

### Always a backup



The EmStat4S is equipped with 500 MB internal storage memory. This means all your measurements can be saved on-board as a backup. All internally stored measurements can be browsed and transferred back to the PC easily using PStace.

Your data is always with your instrument wherever you take it.

## Supported Techniques

The EmStat4S supports the following electrochemical techniques:

### Voltammetric techniques

- |                            |     |
|----------------------------|-----|
| ▪ Linear Sweep Voltammetry | LSV |
| ▪ Cyclic Voltammetry       | CV  |
| ▪ Fast Cyclic Voltammetry  | FCV |
| ▪ AC Voltammetry           | ACV |

### Pulsed techniques

- |                                  |     |
|----------------------------------|-----|
| ▪ Differential Pulse Voltammetry | DPV |
| ▪ Square Wave Voltammetry        | SWV |
| ▪ Normal Pulse Voltammetry       | NPV |

These methods can all be used in their stripping modes which are applied for (ultra-) trace analysis.

### Amperometric techniques

- |                                 |     |
|---------------------------------|-----|
| ▪ Chronoamperometry             | CA  |
| ▪ Zero Resistance Amperometry   | ZRA |
| ▪ Chronocoulometry              | CC  |
| ▪ MultiStep Amperometry         | MA  |
| ▪ Fast Amperometry              | FAM |
| ▪ Pulsed Amperometric Detection | PAD |

### Galvanostatic techniques

- |                              |     |
|------------------------------|-----|
| ▪ Linear Sweep Potentiometry | LSP |
| ▪ Chronopotentiometry        | CP  |
| ▪ MultiStep Potentiometry    | MP  |
| ▪ Open Circuit Potentiometry | OCP |

### Other

- |  |            |
|--|------------|
| ▪ Mixed Mode   | MM         |
| ▪ Potentiostatic and Galvanostatic Impedance spectroscopy at fixed frequency or frequency scan vs <ul style="list-style-type: none"> <li>○ fixed potential or fixed current</li> <li>○ scanning potential or scanning current</li> <li>○ time</li> </ul> | EIS/GEIS   |
| ▪ Fast EIS/GEIS<br>Very low interval fixed-frequency measurements  | FEIS/FGEIS |

MethodSCRIPT™ allows for developing custom techniques. See page 12 for more information.



## Measurement Specifications

The following table shows limits for some technique-specific parameters.

	Parameter	Min	Max
All techniques (unless otherwise specified)	▪ Conditioning time	0	4000 s
	▪ Deposition time	0	4000 s
	▪ Equilibration time	0	4000 s
	▪ Step potential	LR: 0.100 mV HR: 0.183 mV	250 mV
	▪ N data points	3	1 000 000
▪ NPV ▪ DPV	▪ Scan rate	LR: 0.1 mV/s (100 $\mu$ V step) HR: 0.1 mV/s (183 $\mu$ V step)	1 V/s (5 mV step)
	▪ Pulse time	0.4 ms	300 ms
▪ SWV	▪ Frequency	1 Hz	1250 Hz
▪ LSV ▪ CV	▪ Scan rate	LR: 0.01 mV/s (100 $\mu$ V step) HR: 0.01 mV/s (183 $\mu$ V step)	500 V/s (200 mV step)
▪ FCV	▪ Scan rate	LR: 0.1 mV/s (100 $\mu$ V step) HR: 0.01 mV/s (183 $\mu$ V step)	500 V/s (50 mV step)
	▪ N averaged scans	1	65535
	▪ N equil. scans	0	65535
▪ PAD	▪ Interval time	50 ms	4294 s
	▪ Pulse time	1 ms	1 s
	▪ N data points	3	1 000 000 (> 100 days at 10 s interval)
▪ CA ▪ CP ▪ OCP	▪ Interval time	0.4 ms	4294 s
	▪ Run time	1 ms	> year
▪ MM ▪ MA ▪ MP	▪ N cycles	1	20000
	▪ N levels	1	255
	▪ Level switching overhead time	~1 ms (typical)	-
	▪ Interval time	0.4 ms	4294 s
▪ FAM	▪ Interval time	1 $\mu$ s	60 s
	▪ Run time	3 $\mu$ s	34 days (60 s interval) 50 ms (1 $\mu$ s interval)
	▪ N data points	3	50000
▪ Fast EIS	Interval time between points at fixed frequency	~1 ms (typical)	-

## System Specifications

General		
	LR	HR
▪ dc-potential range	±3 V	±6 V
▪ compliance voltage	±5 V	±8 V
▪ maximum current	±30 mA	±200 mA
▪ max. data acquisition rate	1M samples/s	
▪ control loop bandwidth (stability setting)	320 Hz, 3.2 kHz, 30 kHz or 570 kHz	
▪ current follower bandwidth	23 Hz in 1 nA and 10 nA range 2.3 kHz in 100 nA and 1 uA range 230 kHz in 10 uA and 100 uA range > 500 kHz in ranges 1 mA and higher	

Potentiostat (controlled potential mode)		
	LR	HR
▪ applied potential resolution	100 µV	183 µV
▪ applied potential accuracy	≤ 0.2% ±1 mV offset	
▪ current ranges	1 nA to 10 mA (8 ranges)	100 nA to 100 mA (7 ranges)
▪ measured current resolution	0.009% of CR (92 fA on 1 nA range)	
▪ measured current accuracy	< 0.2% of current ±20 pA ±0.2% of range	< 0.2% of current ±0.2% of range

Galvanostat (controlled current mode)		
	LR	HR
▪ current ranges	10 nA, 1 uA, 100 uA, 10 mA (4 ranges)	1 uA, 100 uA, 10 mA, 100 mA (4 ranges)
▪ applied dc-current	±3 * CR (current range)	
▪ applied dc-current resolution	0.01% of CR	0.0183% of CR
▪ applied dc-current accuracy	< 0.4% of current ±20 pA ±0.2% of range	< 0.4% of current ±0.2% of range
▪ potential ranges	50 mV, 100 mV, 200 mV, 500 mV, 1 V	
▪ measured dc-potential resolution	96 µV at ±3 V (1 V range) 48 µV at ±1.5 V (500 mV) 19.2 µV at ±0.6 V (200 mV) 9.6 µV at ±0.3 V (100 mV) 4.8 µV at ±0.150 V (50 mV)	193 µV at ±6 V (1 V range) 96.5 µV at ±3 V (500 mV) 38.5 µV at ±1.2 V (200 mV) 19.3 µV at ±0.6 V (100 mV) 9.65 µV at ±0.3 V (50 mV)
▪ measured dc-potential accuracy	≤ 0.2% potential, ±1 mV offset	

## FRA / EIS (impedance measurements)

▪ frequency range	10 $\mu$ Hz to 200 kHz
▪ ac-amplitude range	1 mV to 900 mV rms, or 2.5 V p-p

## GEIS (galvanostatic impedance measurements)

▪ frequency range	10 $\mu$ Hz to 100 kHz
▪ ac-amplitude range	0.9 * CR A rms

## Electrometer


▪ electrometer amplifier input	> 1 T $\Omega$ // 10 pF
▪ bandwidth	500 kHz

## Other


	LR	HR
▪ electrode connections	WE, RE, CE, and ground, with 2 mm banana plugs	WE, RE, CE, S and ground, with 2 mm banana plugs
▪ power consumption	0.75 W @ 1 mA (WE) 1.25 W @ 30 mA (WE)	1.25 W @ 1 mA (WE) 1.6 W @ 30 mA (WE) 3.75 W @ 200 mA (WE)
▪ power + communication	USB-C	
▪ housing	aluminium body: 7.2 x 5.5 x 2.6 cm	
▪ weight	~130 g	
▪ internal storage space	500 MB, equivalent to >15M datapoints or ~1000 measurement files (whichever comes first)	



## EmStat4S LR EIS Accuracy Contour Plot

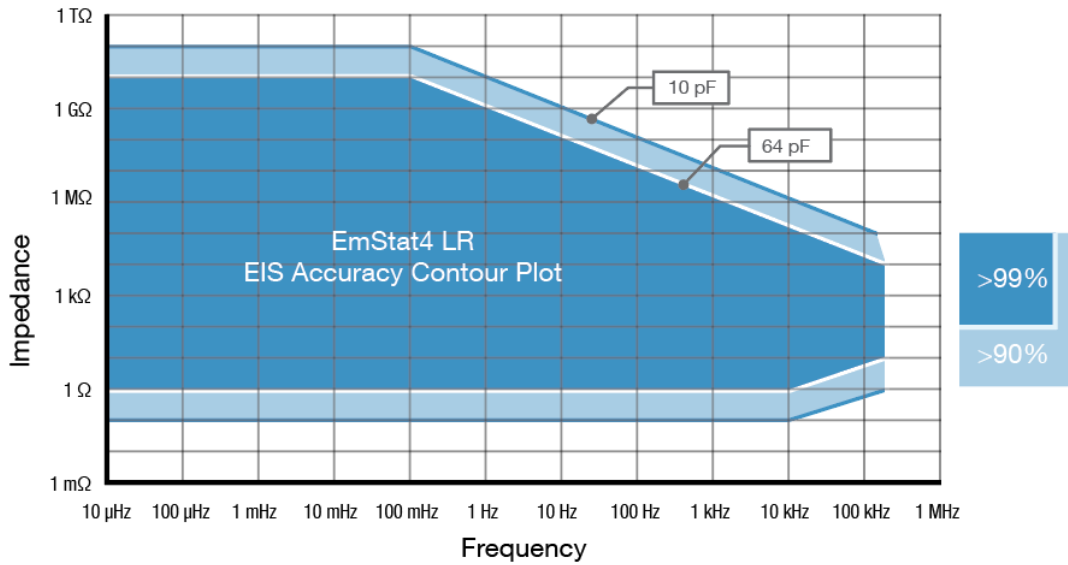

## EmStat4S HR EIS Accuracy Contour Plot

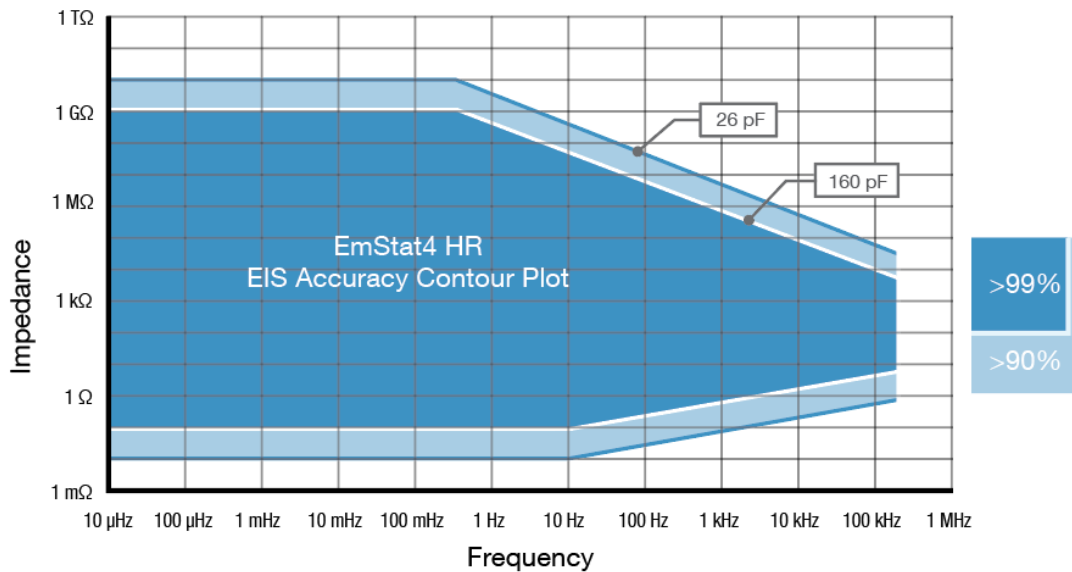

### Note

The accuracy contour plots were determined with an ac-amplitude of  $\leq 10$  mV rms for all limits, except for the high impedance limit, which was determined using an ac-amplitude of 250 mV. The standard 1 meter cell cables were used. Please note that the true limits of an impedance measurement are influenced by all components in the system, e.g. connections, the environment, and the cell.



## Standard EmStat4S Kit

A standard EmStat4S kit includes a rugged carrying case with:

- EmStat4S LR or HR
- USB-C cable
- USB-C splitter cable for extra power (EmStat4S HR only)
- 1 meter cell cable with 2 mm banana pins
- 4 or 5 croc clips
- Dummy Cell

Also included:

- PSTrace software for Windows (on USB drive)
- Manual (hardcopy)
- Quick Start document
- Calibration report



## PSTrace: Software for Windows


PSTrace is designed to get the most out of your instrument right after installation, without going through a long learning period. It has three modes: the Scientific mode which allows you to run all the techniques our instruments have to offer, and two dedicated modes for Corrosion analysis and the Analytical Mode. The Analytical Mode is designed for use with (bio)sensors and allows you to do concentration determinations. Extensive help files and prompts guide the user through a typical analysis.

The screenshot displays the PSTrace software interface. On the left, there is a 'Measurement' panel with settings for 'Cyclic Voltammetry', including current ranges and pretreatment/post-measurement options. The central plot shows 'Current/µA' vs 'Potential/V' with three scans. A 'Selected measurement' dialog box is open, showing details for 'Cyclic Voltammetry' and options to generate new curves. On the right, there is a 'Session data' panel and a 'Find peaks' panel. Three blue callout boxes provide instructions: 'Switch between plots if curves with different units are available.' (pointing to the plot area), 'Setup your measurement easily and get immediate feedback on validity of parameters.' (pointing to the measurement settings), and 'Click on a measurement for detailed information or generating new curves.' (pointing to the measurement dialog). A fourth callout box says 'Quickly toggle the visibility of curves or groups of curves.' (pointing to the session data panel).

## Scripting

The intuitive script editor allows for easily creating a sequence of measurements or other tasks, by means of dragging and dropping actions in a list.

The screenshot shows the 'Script' editor in PSTrace. It features a 'Common' tab with a list of actions such as 'Cell', 'Repeat', 'SetCurrent', 'SetPotential', 'Wait', 'Measurement', and 'FindPeaks'. A script is being built with the following steps: 'Cell', 'Repeat' (1), 'SetPotential' (1.000 V), 'Wait' (5 seconds), 'Measurement' (PSDiffPulse (DPV)), and 'FindPeaks'. The 'FindPeaks' panel on the right is configured with 'Mode: Use window', 'Number of peaks: 1', and a window for Peak 1 from -0.200 to 0.200. The output path is shown as 'C:\Users\Niels van Velzen\CloudStation\PSData\scriptout...'. A 'Remove' button and navigation arrows are at the bottom.

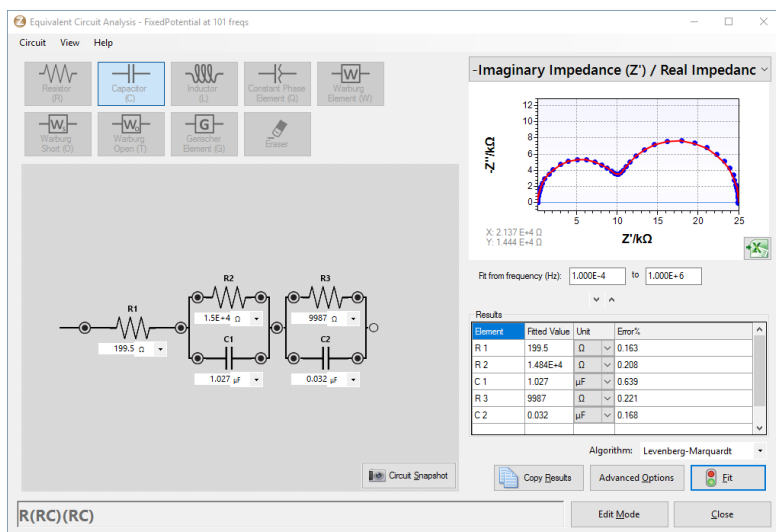


Use the graphical editor to draw the equivalent circuit or enter the CDC directly.

Corrosion mode for Tafel plot analysis and other corrosion data analysis.

## Other functions in PSTrace

- Concentration determination
- Advanced peak search algorithms
- Open your data in Origin and Excel with one click of a button
- Save all available curves, measurement data and methods to a single file
- Load measurements from the internal storage
- Direct validation of method parameters
- Run custom MethodSCRIPTS™

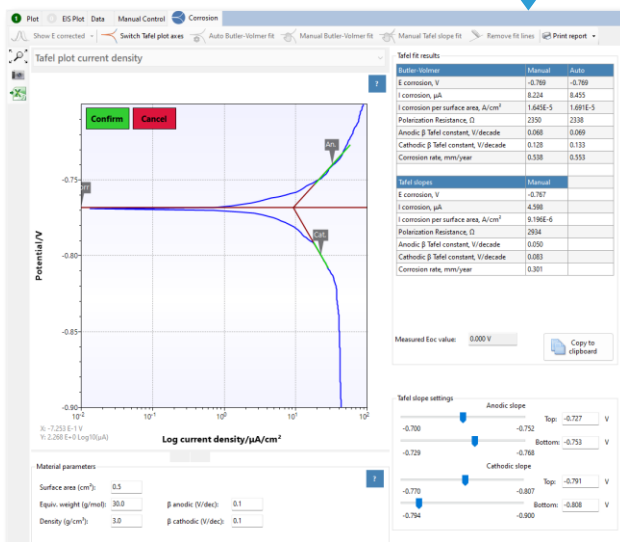

## Integration with third party software

- Excel
- Origin
- Matlab
- ZView



## Minimum System Requirements

- Windows 7, 8, 10 or 11
- 1 GHz or faster 32-bit (x86) or 64-bit (x64) processor
- 2 GB RAM (32-bit) or 4 GB RAM (64-bit)
- Screen resolution of 1280 x 800 pixels

➤ See for more information:  
[www.palmsens.com/pstrace](http://www.palmsens.com/pstrace)

## EmStat4S works with MethodSCRIPT™

The MethodSCRIPT™ scripting language is designed to integrate our instruments and potentiostat (modules) effortlessly in your hardware setup, product, or experiment.

MethodSCRIPT™ gives you full control over your potentiostat. The simple script language is parsed on-board the instrument and allows for running all supported electrochemical techniques, making it easy to combine different measurements and other tasks.


MethodSCRIPT can be generated, edited, and executed in PStRace.

MethodSCRIPT features include:


- Use of variables
- (Nested) loops and conditional logic support
- User code during a measurement iteration
- Exact timing control
- Simple math operations on variables (add, sub, mul, div)
- Digital I/O, for example for waiting for an external trigger
- Logging results to internal storage or external SD card
- Reading auxiliary values like pH or temperature
- and many more...

```
1 e
2 var c
3 var p
4 #Select bandwidth of 40 for 10 points per second
5 set_max_bandwidth 40
6 #Set current range to 1 mA
7 set_range ba 1m
8 #Enable autoranging, between current of 100 uA and 1 mA
9 set_autoranging ba 100u 1m
10 #Turn cell on for measurements
11 cell_on
12 #equilibrate at -0.5 V for 5 seconds, using a CA measurement
13 meas_loop_ca p c -500m 500m 5
14 pck_start
15 pck_add p
16 pck_add c
17 pck_end
18 endloop
19 #Start LSV measurement from -0.5 V to 1.5 V, with steps of 10 mV
20 #and a scan rate of 100 mV/s
21 meas_loop_lsv p c -500m 1500m 10m 100m
22 #Send package containing set potential and measured WE current.
23 pck_start
24 pck_add p
25 pck_add c
26 pck_end
27 #Abort if current exceeds 1200 uA
28 if c > 1200u
29 abort
30 endloop
31 #Turn off cell when done or aborted
32 on_finished:
33 cell_off
34
```







[Online support on MethodSCRIPT](#)

 Write your own software and integrate (generated) MethodSCRIPTs. No libraries needed.

MethodSCRIPT is parsed on-board the instrument. No DLLs or other type of code libraries are required for using MethodSCRIPT™

  
MethodSCRIPT™

Code examples are available for:

-  C/C++
-  ARDUINO
-  Swift
-  Xamarin
-  python™
-  Java

➤ See for more information:  
[www.palmsens.com/methodscript](http://www.palmsens.com/methodscript)

## Software Development Kits for .NET

Develop your own application in no time for use with any PalmSens instrument or potentiostat (module). Our SDKs are free of charge.



There are three PalmSens Software Development Kits (SDKs) for .NET. Each SDK can be used with any of our instruments or OEM potentiostat modules to develop your own software. The SDK's come with a set of examples that shows how to use the libraries. PalmSens SDKs with examples are available for the following .NET Frameworks:

- WinForms
- Xamarin (Android)
- WPF

Each SDK comes with code examples for:

- Connecting
- Running measurements and plotting data
- Manual control of the cell
- Accessing and processing measured data
- Analyzing and manipulating data
- Peak detection
- Equivalent Circuit Fitting on impedance data
- Saving and loading files

```
/// <summary>
/// Initializes the EIS method.
/// </summary>
/// </summary>
1 reference
private void InitMethod()
{
    _methodEIS = new ImpedimetricMethod();
    _methodEIS.ScanType = ImpedimetricMethod.enumScanType;
    _methodEIS.Potential = 0.0f; //0.0V DC potential
    _methodEIS.Eac = 0.01f; //0.01V RMS AC potential at
    _methodEIS.FreqType = ImpedimetricMethod.enumFrequency;
    _methodEIS.MaxFrequency = 1e5f; //Max frequency is
    _methodEIS.MinFrequency = 10f; //Min frequency is
    _methodEIS.nFrequencies = 11; //Sample at 11 different
    _methodEIS.EquilibrationTime = 1f; //Equilibrates
    _methodEIS.Ranging.StartCurrentRange = new CurrentRange;
    _methodEIS.Ranging.MinimumCurrentRange = new CurrentRange;
    _methodEIS.Ranging.MaximumCurrentRange = new CurrentRange;
}
```

➤ See for more information:  
[www.palmsens.com/sdk](http://www.palmsens.com/sdk)

## PStouch: App for Android



PStouch is an app for Android devices compatible with all PalmSens, EmStat and Sensit potentiostats. Your smartphone or tablet connects with the EmStat4S via USB. This requires your device to support USB OTG.

### PStouch features:

- Setting up and running measurements
- Loading and saving measured curves
- Analyzing and manipulating peaks
- Sharing measurement data directly via any service like email or Dropbox
- Concentration determination by means of Standard Addition or Calibration Curve
- Support for PalmSens accessories such as a Multiplexer or Stirrer
- All method and curve files are fully compatible with PStace software for Windows.

➤ See for more information:  
[www.palmsens.com/pstouch](http://www.palmsens.com/pstouch)

Please do not hesitate to contact PalmSens for more details: [info@palsens.com](mailto:info@palsens.com)

**PalmSens BV**  
**The Netherlands**  
[www.palsens.com](http://www.palsens.com)

**DISCLAIMER**

Changes in specifications and typing errors reserved. Every effort has been made to ensure the accuracy of this document. However, no rights can be claimed by the contents of this document.