


OG Carbon New Silver


Screen Printed Carbon Electrodes

Droplet placed here

Plug in here

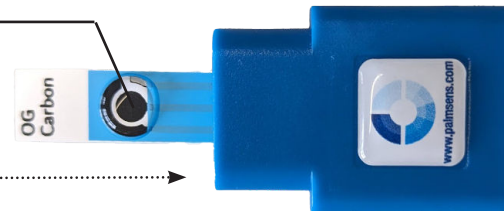



Disposable and high quality screen printed carbon electrodes with good reproducibility. Perfect for biosensor research and development.

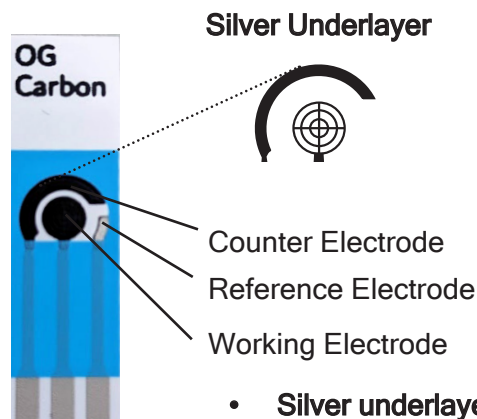


- **Silver underlayer** decreases track resistance
- Carbon (graphite) counter and working electrode
- Silver/silver chloride reference electrode
- 4 mm working electrode diameter
- Flexible PET support material
- 10 mm x 33 mm x 0.25 mm
- Compatible with a range of potentiostats and connectors

Available in rows of 25

Easily cuttable

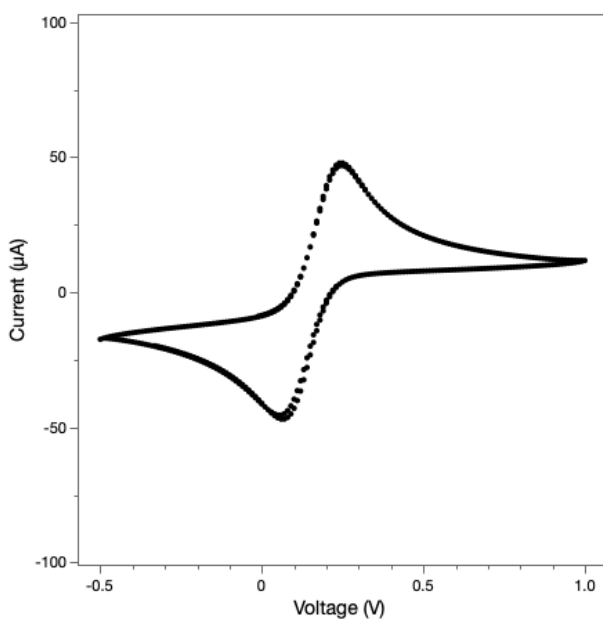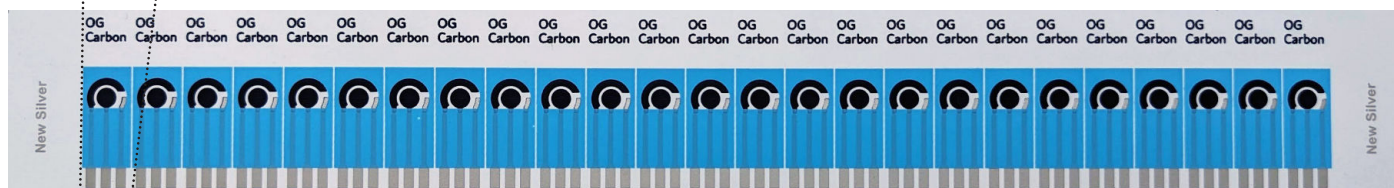

Cyclic voltammetry of 3 OG Carbon New Silver electrodes in $\text{Fe}(\text{CN})_6^{3-/4-}$

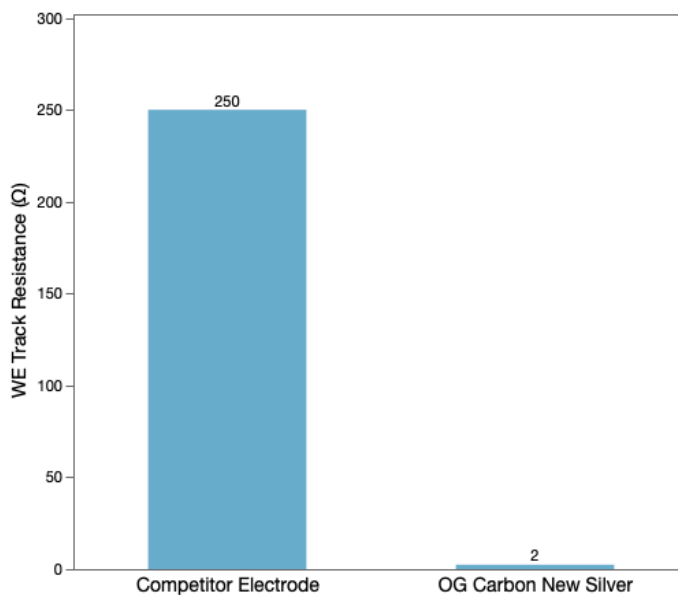

Working track resistance for OG Carbon new silver electrodes vs competitor electrodes

More information available at ogcarbon.co.uk